


PM-COMFORT™

MATTRESS



- HIGH COMFORT MATTRESS
- OPTIMIZED FOR XRAY
- HYGIENIC AND DURABLE

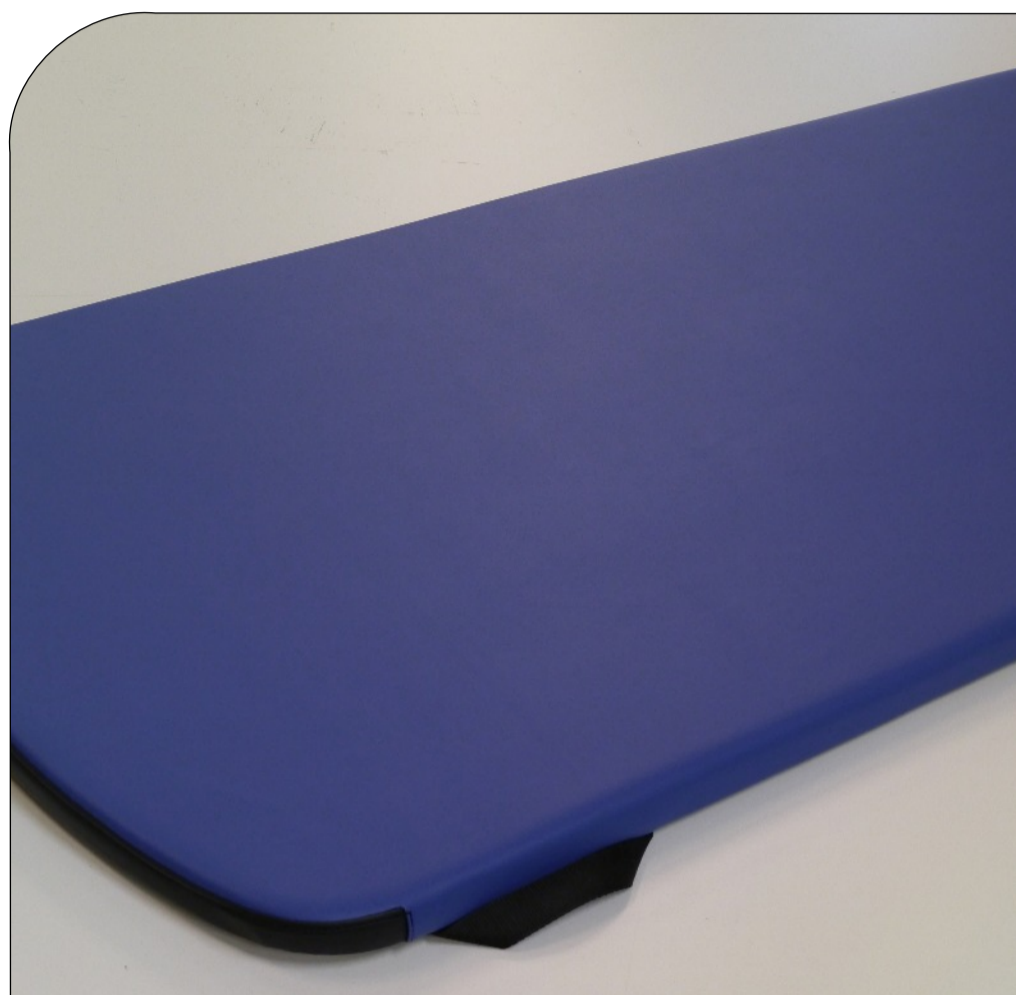

PM Comfort™ is a pressure relieving mattress for patient including intervention / angio and convetional labs. The mattress can be equipped with handles and lowfriction underside, which facilitates positioning. Delivered to the leading radiology companies. To the right HighPos™ Shaped - wedge cushion for leg.



ESON COMFORT
Idé-Möbler i Landeryd AB
Centrumvägen 7
SE-314 94 LANDERYD

Tel +46 (0)371-40051
info@esoncomfort.com

Copyright 2016 Eson Comfort / Idé-Möbler i Landeryd AB.

All rights reserved

PRODUCTINFORMATION

The mattress consists of a 25mm polyurethane foam with a density of approximately 52kg / m³. The core consists of a 9mm's EVA plate with a volume weight of about 50kg / m³. The cover is made of syntetic leather and can handle cigarette and ignition test according to Swedish a minimum standard, without the addition of flame retardants. Production without CFCs. PM Comfort™ has a technological function warranty of 1 year.

ART.NO
PMCOMF
PMCOMFH

NAME
Patient mattress Width x Length x 3cm
Patient mattress Width x Length x 3 cm with Handles

NOTE! Drawing or template must always be enclosed the order of molded mattresses.

Option / spare parts
Special shape or special dimensions.

Washing
The cover washed with spirit of soap solution.

CE Registrationcode 344-2757-35331

ESON COMFORT
WWW.ESONCOMFORT.COM
"WE INVENT THE FUTURE"